


# 2200G SERIES

## N95 HEALTHCARE AND RESPIRATORY SURGICAL MASK



Available in 4 color coded sizes

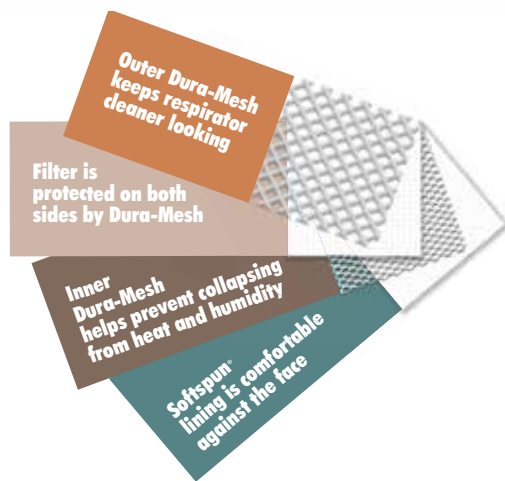

2210GN95 XS



2211GN95 S



2212GN95 M/L



2217GN95 LP

### A FIT FOR EVERYONE

Available in 4 color-coded (strap) sizes for today's diverse workforce.

### NON-LATEX RESPIRATOR AND STRAPS

The non-latex straps don't irritate or require prestretching, plus they are securely attached.

### EASY FIT, MORE COMFORT

The molded nose bridge with soft foam nose cushion provides the wearer an easy comfortable fit. There's no metal noseband to adjust and the seal is maintained without pressure points.

### DURA-MESH® SHELL LASTS LONGER

The exclusive Dura-Mesh shell protects the filter media so it stays cleaner-looking longer. The shell also resists collapsing in hot, humid environments so the respirator holds its shape longer. The result is fewer respirators used and lower costs.

### FEATURES

- Meets CDC guidelines for protection against TB.
- NIOSH certified to have a filter efficiency of 95% or greater against non-oil based particulates.
- Bacterial filter efficiency >99.9%.
- Non-Latex respirator and straps.
- Dura-Mesh shell resists collapsing in regular use and even under hot and humid conditions.
- Facepiece is naturally contoured without pressure points.
- Soft foam nose cushion for added comfort.
- Softspun® inner lining is soft to the touch.
- Can be used in emergency situations by paramedics.
- Suggested for stockpiling.
- Meets heat and flame resistance in accordance with ANSI/ISEA 110-2009 Section 7.11.1.
- 100% PVC-Free for a greener alternative.

# 2200G SERIES

## N95 HEALTHCARE AND RESPIRATORY SURGICAL MASK

Product Code	Size	Valve Type	Strap Material	Display Quantity	Case Quantity
2210GN95-XS	XS	None	Non-Latex	20/Box	12 Boxes
2211GN95-S	S	None	Non-Latex	20/Box	12 Boxes
2212GN95-M/L	M/L	None	Non-Latex	20/Box	12 Boxes
2217GN95-LP	LP	None	Non-Latex	20/Box	12 Boxes

XS = Extra Small S = Small M/L = Medium/Large LP = Low Profile Nose Bridge

### WARNING TO USER

1. Non-latex straps are used in the headbands of these respirators. However, individuals highly sensitive to natural rubber latex may have an allergic reaction.
2. This product does not eliminate the risk of contracting any disease or infection. Change immediately if the respirator comes in contact with blood or fluids, including body fluids. OSHA has not set a permissible exposure level for airborne biohazards.
3. This product has not been sold with warning or use instructions for personnel involved in industrial or related situations. If you are considering such uses, first contact the Moldex Technical Services Department at +1 (310) 837-6500 or +1 (800) 421-0668, ext. 512/550.

### DESCRIPTION

1. The Moldex Healthcare N95 Particulate Respirator and Surgical Mask is designed to help provide respiratory protection for the wearer. This product has been tested<sup>1</sup> and certified by NIOSH as an N95 respirator and as having a filter efficiency level of 95% or greater against particulate aerosols free of oil. It is fluid resistant,<sup>2</sup> disposable and may be worn in surgery or throughout the hospital.
2. This product contains no components made from natural rubber latex.

### INTENDED USE

The various models of Moldex Type N95 Healthcare Particulate Respirators and Surgical Masks meet CDC Guidelines for TB Exposure Control within healthcare facilities. These devices are also intended to be worn by healthcare personnel during surgical procedures to protect both the patient and healthcare personnel from the transfer of microorganisms, body fluids and particulate material.

### CONTRAINDICATIONS

1. Not for use with beards, or other facial hair that prevents direct contact between the face and sealing surface of the respirator.
2. Eyewear must not prevent direct contact between the face and sealing surface of the respirator.
3. Not to be used on children.


<sup>1</sup> Tested in accordance with NIOSH 42 CFR 84

<sup>2</sup> Passed ASTM F 1862 @ 160mm Hg




**MOLDEX-METRIC, INC.**  
 10111 Jefferson Blvd.  
 Culver City, CA 90232  
 TEL: +1 (800) 421-0668 or  
 +1 (310) 837-6500  
 FAX: +1 (310) 837-9563  
 E-mail: sales@moldex.com  
 www.moldex.com

**CANADIAN CUSTOMER SERVICE**  
 Tel: +1 (800) 421-0668, Ext. 550/512  
 Fax: +1 (310) 837-9563



Moldex-Metric, Inc.  
 Culver City, California 90232  
 +1 (800) 421-0668



THESE RESPIRATORS ARE APPROVED ONLY IN THE FOLLOWING CONFIGURATIONS:

TC	PROTECTION <sup>1</sup>	RESPIRATOR				CAUTIONS AND LIMITATIONS <sup>2</sup>
		2210GN95-XS	2211GN95-S	2212GN95-M/L	2217GN95-LP	
84A-4339	N95	X	X	X	X	A B C J M N O P

1. Protection
 

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply
2. Cautions and Limitations
  - A- Not for use in atmospheres containing less than 19.5 percent oxygen.
  - B- Not for use in atmospheres immediately dangerous to life or health.
  - C- Do not exceed maximum use concentrations established by regulatory standards.
  - J- Failure to properly use and maintain this product could result in injury or death.
  - M- All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
  - N- Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
  - O- Refer to User's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
  - P- NIOSH does not evaluate respirators for use as surgical masks.

### RESTRICTIONS

1. Before use, a written respiratory protection program in accordance with 29 CFR 1910.134 must be implemented.
2. This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.
3. Prior to each use, carefully inspect the entire respirator, including filter media and strap attachment area for tears and damage. Staple or strap attachment perforations do not effect NIOSH approval.
4. If respirator is damaged, distorted, a proper fit cannot be obtained or breathing becomes difficult, leave contaminated area as soon as possible and replace respirator.
5. If the respirator comes in contact with blood or fluids, including body fluids, leave contaminated area as soon as possible and discard and replace the respirator.
6. Do not alter, modify or abuse this respirator.
7. Store unused respirators in box/bag in a clean, dry, non-contaminated area.
8. Dispose of respirator according to your employer's policy and local regulations.
9. Do not reuse or store for reuse or hang around neck unless your employer specifically authorizes reuse.
10. When used for surgical procedures, discard after every use.
11. Use respirator before the "use by" expiration date printed on box/bag.

### LIMITED WARRANTY IMPORTANT NOTICE TO PURCHASER

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Moldex will be met by products sold. Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Moldex, products proved to be defectively manufactured. Except as provided above, Moldex shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

### MOLDEX TECHNICAL ASSISTANCE

For more information call +1 (800) 421-0668 or +1 (310) 837-6500 ext. 512/550 or E-mail: sales@moldex.com, tech@moldex.com or visit www.moldex.com.

### DISTRIBUTED BY:



Printed on Recycled Paper